
**CiD**  
Carbostent & Implantable Devices

INTERVENTIONAL CARDIOLOGY

**flyid's**  
PTCA RAPID EXCHANGE BALLOON CATHETER

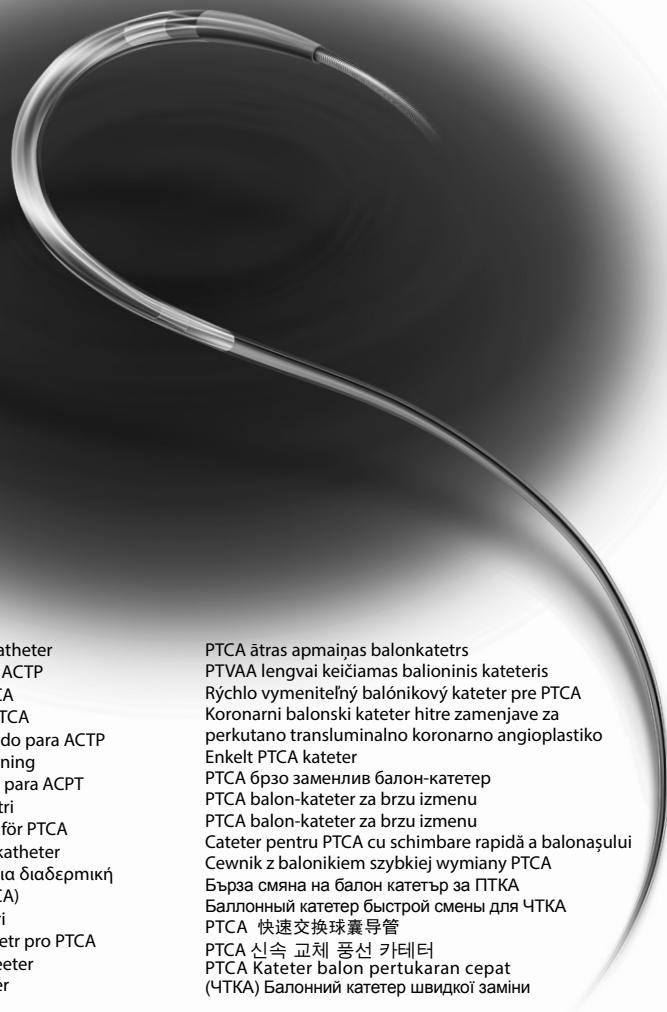PTCA Rapid exchange balloon catheter  
Cathéter à échange rapide pour ACTP  
Schnellwechsel-Katheter für PTCA  
Catetere a scambio rapido per PTCA  
Catéter de balón de cambio rápido para ACTP  
Kateter til PTCA til hurtig udskiftning  
Cateter de balão de troca rápida para ACPT  
Pikavaihettava PCTA-pallokatetri  
Snabbt utbytbar ballongkater för PCTA  
Snel verwisselbare PTCA-ballonkatheter  
Καθετήρας ταχείας εναλλαγής για διαδερμική οτεφανιάσια αγγειορλαστική (PTCA)  
PTCA hızlı değişim balon kateteri  
Rychle měnitelný balónkový katetr pro PTCA  
PTKA kiirvahetusballooniga kateeter  
PTCA gyors csereballonos katéter

PTCA ātras apmaiņas balonkatetrs  
PTVAA lengvai keičiamas balioninis kateteris  
Rýchlo vymenitelný balónkový kateter pre PTCA  
Koronarni balonski kateter hitre zamenjave za perkutano transluminalno koronarno angioplastiko  
Enkelt PTCA kateter  
PTCA брзо заменлив балон-катетер  
PTCA balon-kateter za brzu izmenu  
PTCA balon-kateter za brzu izmenu  
Cateter pentru PTCA cu schimbare rapidă a balonașului  
Cewnik z balonkiem szybkiej wymiany PTCA  
Бърза смяна на балон катетър за ПТКА  
Баллонный катетер быстрой смены для ЧТКА  
PTCA 快速交换球囊导管  
PTCA 신속 교체 풍선 카테터  
PTCA Kateter balon pertukaran cepat (ЧТКА) Балонний катетер швидкої заміни

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## ENGLISH

### DESCRIPTION

The FLUYDO PEGASO percutaneous transluminal coronary angioplasty (PTCA) catheter (shown in the figure) is a rapid exchange type catheter. Balloon diameters at different pressures (Compliance) are shown in the table 2.

The distal portion of the catheter has a hydrophilic coating and consists of two lumens: one for balloon inflation and deflation, the other for guidewire advancement and retraction.

Two radiopaque markers located at both ends of the balloon permit accurate placement of the balloon across the stenosis under fluoroscopy (NOTE: In models with a nominal diameter of 1.5 mm, the balloon has a single radiopaque marker placed in the centre).

The proximal portion of the catheter, a stainless steel hypotube, contains the lumen for balloon inflation and deflation.

Two depth markers, one at 90 and the other at 100 cm from the distal tip, aid in determining when the balloon is exiting the guiding catheter, in the cases of brachial or femoral/radial approach, respectively.

The proximal end of the catheter has a female Luer lock for attachment to an inflation device.

The Manufacturer produces the FLUYDO PEGASO directly and carries out the quality control, both during the manufacturing process and on the finished product, in accordance with the appropriate Good Manufacturing Standards.

The FLUYDO PEGASO catheter is supplied sterile, non-pyrogenic and individually packed in a pouch that must not be placed in a sterile field.

The sterilization made by the manufacturer uses an ethylene oxide and CO<sub>2</sub> mixture.

Sterility is guaranteed while the package remains intact and until the expiry date printed on the packaging (EXPIRY DATE).

### Contents:

- One Fluydo Pegaso balloon dilatation catheter
- One compliance table

### Storage:

Store in a cool, dry place, away from sunlight.

### INTENDED USE

The FLUYDO PEGASO catheter is indicated for balloon dilatation of the stenotic portion of a coronary artery or bypass graft stenosis for the purpose of improving myocardial perfusion.

### INDICATIONS

The FLUYDO PEGASO catheter is indicated in the treatment of obstructive lesions of native coronary arteries or aortocoronary bypasses in patients with the requirements for percutaneous transluminal coronary angioplasty (PTCA) and for post dilation of coronary stents.

### CONTRAINDICATIONS

The FLUYDO PEGASO catheter is contraindicated in patients with a previous diagnosis of coronary artery spasm in the absence of a significant stenosis. The device is not intended to be used in an unprotected left main coronary artery.

### WARNINGS

- The device is intended for single use only. Do not reuse, reprocess or resterilise. This could compromise device performance and cause contamination risks of the device and/or infections to the patient, inflammation and patient-to-patient transmission of infective diseases.
- Do not use a catheter if its packaging has been opened or damaged, or if the "Use by" date has expired.
- PTCA in patients thought to be unsuitable candidates for coronary bypass requires very careful consideration and possible haemodynamic support during PTCA, since treatment of this type of patient involves particular risk.
- The device must be handled in such a way as to avoid contact with metal or abrasive instruments, since this could damage the catheter and compromise its performance.
- The catheter must not be used if the shaft of the catheter has kinks or bends; in this case, do not attempt to straighten the catheter.
- The catheter must be handled under fluoroscopic control using X-ray equipment that produces high-quality images.
- The PTCA catheter should only be advanced over a guidewire. Do not remove the guidewire from the catheter during use.

- Do not advance or retract the catheter unless the balloon is fully deflated and under vacuum.
- If, during any stage of the procedure, unusual resistance occurs, do not force the catheter: withdraw the guiding catheter and the FLUYDO PEGASO catheter together, as if they were one unit. Application of excessive force and/or incorrect handling could damage the catheter.
- Do not pressurise the catheter until the balloon has been positioned at the site of the lesion.
- Never use air or other gaseous substances to inflate the balloon.
- Use of an inflation device with manometer is recommended. While inflating the device do not exceed the nominal value of Rated Burst Pressure.
- The balloon diameter must never exceed the diameter of the coronary artery just proximal and distal to the stenosis. Vascular damage could occur.
- A surgery team should be available for a possible intervention.
- Do not attempt to reclean or resterilise devices that have been in contact with blood or body tissues. Used devices are to be disposed of as dangerous medical waste with risk of infection.
- When using accessories (guide catheter, guidewire, haemostatic valve), follow the Manufacturer's instructions.

### PRECAUTIONS

- Before the angioplasty procedure, examine the catheter to check that it operates correctly (no kinks, bends or other damage) and ensure that the dimensions are correct for the specific procedure in which it is to be used.
- The catheter must only be used by, or under the supervision of, medical personnel suitably qualified in performing percutaneous transluminal coronary angioplasty procedures.
- Before inserting the catheter, administer the appropriate anticoagulant and coronary vasodilatory treatment.
- In the case of simultaneous use of two balloon catheters in one guide catheter, care should be taken when introducing and removing guide wires and balloon catheters to avoid entanglement.

### POTENTIAL ADVERSE EVENTS

Possible adverse effects of PTCA include, but are not limited to:

- Myocardial infarction
- Cardiac arrhythmia, including ventricular fibrillation
- Unstable angina
- Dissection, perforation, rupture of the artery
- Embolism
- Restenosis of the dilated vessel
- Allergic reactions (to contrast medium and drugs used during the procedure)
- Vessel occlusion
- Coronary spasm
- Hypo/hypertension
- Infection
- Haemorrhage or haematoma at the access site
- Arteriovenous Fistula
- Pseudoaneurysm (at site of catheter insertion)
- Thrombosis.

These complications may also require emergency coronary bypass (CABG) surgery or emergency medical treatment; these treatments may have fatal results

### DIRECTIONS FOR USE

#### Auxiliary material (not included in the box)

- Guiding catheter(s) with a 5F diameter (inner diameter 1.47 mm) or larger
- Syringes
- Guidewires 0.014 inches in diameter (0.356 mm) or less
- Inflation device with manometer
- Introducer

#### Catheter preparation

- a) The inflation diameter of the balloon catheter must not exceed the diameter of the coronary artery just proximal and distal to the stenosis.
- b) After having checked that the pack has not been damaged, remove the catheter from it in sterile conditions, without bending it, and place it on a sterile surface.
- c) Check that the catheter has no kinks, bends or other damage. Remove the distal stiffening wire and the balloon protection.
- d) Wash the guidewire lumen with heparinised solution.
- e) Purge any air from the catheter as follows:
  1. Fill a 20-30 cc syringe with 3-5 ml contrast mixture (50% contrast medium and 50% sterile saline solution) and purge the remaining air.

2. Attach the syringe to the Luer connection.
  3. With the tip of the catheter pointing downwards, aspirate for 30 seconds. Slowly release the plunger of the syringe and allow the contrast medium to flow into the distal portion of the catheter.
  4. Disconnect the syringe from the luer connection and remove all air from the syringe.
  5. Repeat step 3, aspirate the device for 10-15 seconds until bubbles no longer appear. Remove the syringe.
- f) Prepare the inflation device according to the manufacturer's instructions. Attach the inflation device to the inflation port to make fluid-to-fluid contact.
- g) Before inserting the catheter, apply negative pressure to the balloon with the inflation device. The negative pressure will provide the lowest possible balloon profile to ease catheter insertion.

#### Inserting the catheter

Guiding catheters with a diameter of 5F (inner diameter 1.47 mm) or more may be used with FLUYDO PEGASO catheters

Coronary guidewires with a diameter of 0.014 inches (0.356 mm) or less may be used with FLUYDO PEGASO catheters. The doctor may choose the stiffness and type of tip according to his own clinical experience.

- a) Clean the exposed portion of guidewire with a mixture of saline solution and heparin, in order to remove traces of blood and contrast liquid.
- b) Ensure that the haemostatic valve is fully open before inserting the catheter.
- c) Thread the distal tip of the catheter over the proximal end of the guidewire until the guidewire exits the catheter body at the guidewire exit port.
- d) Advance the dilatation catheter until the appropriate depth marker aligns with the haemostatic valve hub. This indicates that the balloon is exiting the guiding catheter.
- e) Under fluoroscopy advance the catheter over the guidewire until the catheter radiopaque markers show the balloon is centered in the area to be dilated. Confirm proper positioning of the balloon using fluoroscopy.

**CAUTION: The catheter must not be used if, during handling or insertion, the shaft of the catheter has kinks or bends; in this case, do not attempt to straighten the catheter.**

- f) Kissing balloon: bench tests have shown that two Fluydo Pegaso catheters can be introduced in 6F (or larger) guide catheter.
- g) With the catheter positioned within the stenosis, inject the diluted contrast medium into the inflation lumen of the catheter to inflate the balloon by means of an inflation device. Maintain negative pressure on the balloon between inflations.
- h) The success of the dilation procedure can be determined by injecting contrast medium through the guiding catheter. Ensure the balloon is fully deflated.
- i) Ensure the balloon is fully deflated and remove the PTCA catheter. A device with larger and longer balloon requires a longer deflation time.

**NOTE: Before withdrawal of the PTCA catheter, wipe the guidewire with saline soaked gauze to remove any contrast medium.**

#### **LIABILITY AND WARRANTY**

The Manufacturer guarantees that this device has been designed, manufactured and packaged with the greatest care, using the most appropriate procedures the current state of technology allows. The safety standards integrated into the design and manufacture of the product guarantee its safe use under the aforesaid conditions and for its intended purposes in observation of the precautions listed above. These safety standards are meant to reduce as far as possible, but not eliminate completely, the risks related to the use of the product.

The product must be used only under a specialist physician's care, while taking into account any risks or side effects and complications that might arise from its intended use, as mentioned in other sections of this instructions booklet.

Given the technical complexity, the critical nature of treatment choices and the methods used to apply the device, the Manufacturer cannot be held responsible, either explicitly or implicitly, for the quality of end results subsequent to the use of the device or its effectiveness in resolving a patient's state of infirmity. The end results, in terms of both the clinical status of the patient and the functionality and the lifetime of the device, depend on many factors outside the manufacturer's control, among which are the patient's conditions, the surgical procedure of implantation and application, and the handling of the device after being removed from its package.

In light of these factors, therefore, the Manufacturer is solely responsible for the replacement of any device which, upon delivery, is found to

have manufacturing defects. The buyer shall return the device to the Manufacturer, which reserves the right to inspect the returned device and, at its sole discretion, to determine whether the device is truly defective in manufacturing or material. The warranty consists exclusively in the replacement of the defective device with a Manufacturer device of the same fabrication or equivalent.

The warranty applies under the condition that the device is returned correctly packed to the Manufacturer, and that it is accompanied by a written, detailed report describing the defects claimed, and, if the device has been implanted, stating the reasons for its removal from the patient. On replacement of the device, the Manufacturer shall reimburse the buyer the expenses incurred for replacement of the defective device.

The Manufacturer declines any responsibility for cases of negligence to observe the methods of use and the precautions stated in this instructions booklet and for cases of use of the device after the use before date printed on the package.

Furthermore, the Manufacturer declines any responsibility related to the consequences arising from the treatment choices and the methods of use or the application of the device; the Manufacturer shall, therefore, not be liable for any damage of any nature, material, biological or moral, subsequent to the application of the device.

The agents and the representatives of the Manufacturer are not authorized to modify any of the conditions of this warranty or to assume any further obligations or to offer any guarantees related to this product beyond the terms stated above.

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